



Kwiktrak Configurations



- From the worlds leading ceiling hoist supplier, a patented quick and easy ceiling hoist solution

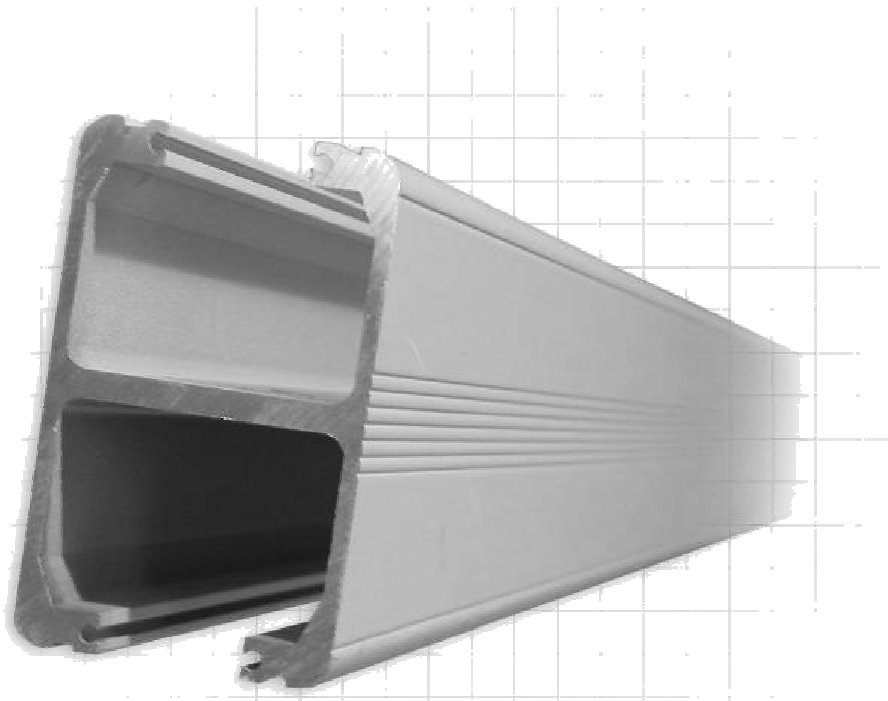

**Kwiktrak
Installations**

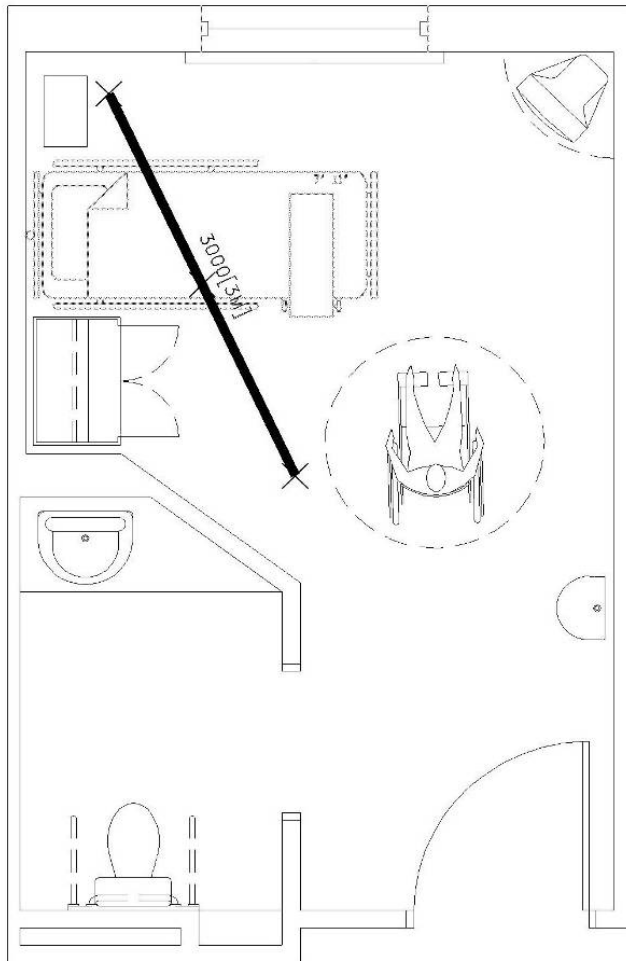

LIST OF MATERIAL		
QTY	PIECE	DESCRIPTION
1	200.12810	3 METERS KWIKTRAK
2	700.11000	END STOPPER
2	200.11060	END CAP
3	700.11100	TRACK BRACKET
1	9100001	V4 - 4 WAYS


Private Room Straight Track:


Configuration Advantages:

- Cost effective use of funds
- Quickest install time
- 4-way lift allows powered movement left and right along the rail
- Return to charge to charger that can be positioned anywhere on rail
- Use of lift for repositioning and turning in bed
- Installed Price £1250

Approved: _____

QTY: 1

 BHM MEDICAL INC.		DISTRIBUTOR: BHM MEDICAL INC.	
REFERENCE #: B040213		1	
REVISION: ---		REP.: CHRISTIAN DEGRACE	
DRAWN BY: F. GUERTIN		DATE (D/M/Y): 13-02-2004	
SCALE: 1=50		CUSTOMER: ---	

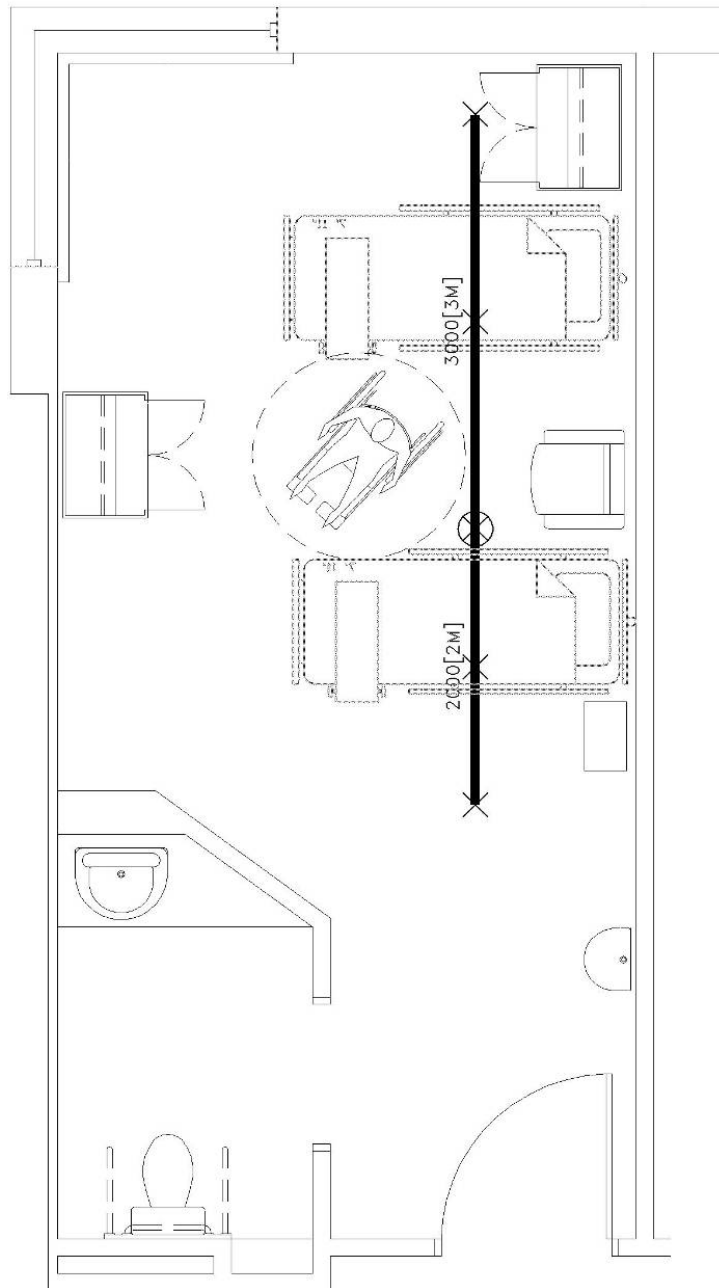

LIST OF MATERIAL		
QTY	PIECE	DESCRIPTION
1	200.12800	2 METERS KWIKTRAK
1	200.12810	3 METERS KWIKTRAK
2	700.11000	END STOPPER
2	200.11060	END CAP
5	700.11100	TRACK BRACKET
2	200.11105	CAM FOR JOINT
2	000.00810	SPRING PIN
1	9100001	V4 - 4 WAYS


Semi-Private Room Straight Track:

Configuration Advantages:

- One lift covers both beds
- Multiple transfer points
- Cost effective use of funds
- 4-way lift allows powered movement left and right along rail
- Return to charge to anywhere on rail
- Use of lift for repositioning and turning in bed
- Installed price £1400

Approved: _____

QTY: 1


 BHM MEDICAL INC.		DISTRIBUTOR: BHM MEDICAL INC.	
REFERENCE #: B040213		3	
REVISION: ---		CUSTOMER: ---	
DRAWN BY: F. GUERTIN		SCALE: 1=50	
		DATE (MM/YY): 13-02-2004	



Straight track installed in Semi-Private room with concrete ceiling. This is accomplished using powder coated exposed uni-strut.



**Kwiktrak
Installations**

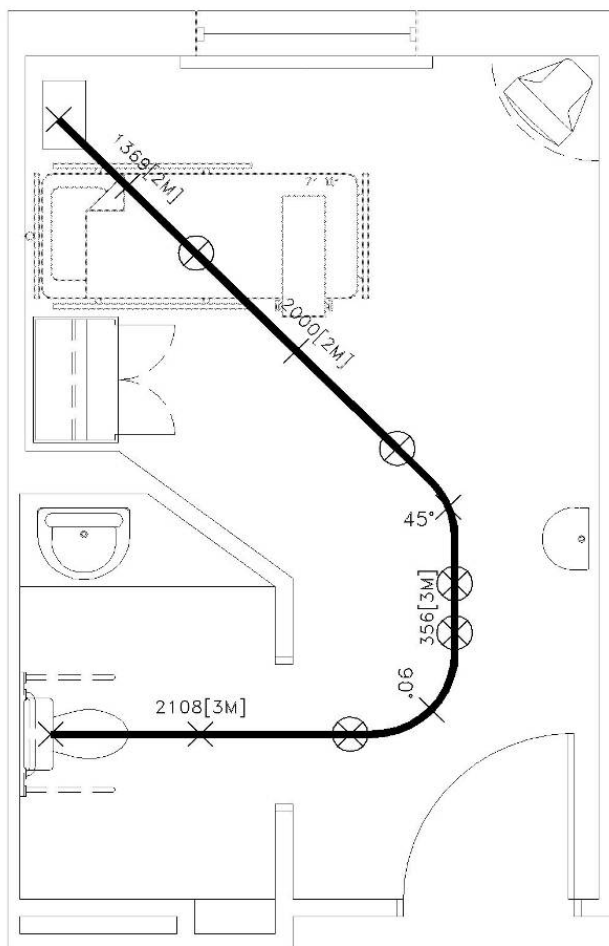

LIST OF MATERIAL		
QTY	PIECE	DESCRIPTION
2	200.12800	2 METERS KWIKTRAK
1	200.12810	3 METERS KWIKTRAK
1	700.12840	45° CURVE KWIKTRAK
1	700.12860	90° CURVE KWIKTRAK
2	700.11000	END STOPPER
2	200.11060	END CAP
12	700.11100	TRACK BRACKET
10	200.11105	CAM FOR JOINT
10	000.00810	SPRING PIN
1	9100001	V4 - 4 WAYS


Private Room Curved Track:

Configuration Advantages:

- Access to washroom with lift
- Multiple transfer points
- 4-way lift allows powered movement left and right along rail
- Use of lift for repositioning and turning in bed
- Installed price guide – approx £1800-£2000

Approved: _____

QTY: 1


 BHM MEDICAL INC.		DISTRIBUTOR:	
REFERENCE #:		B040213	
REVISION:		---	
DRAWN BY:		F. GUERTIN	
SCALE:		1=50	
DISTRIBUTOR:		BHM MEDICAL INC.	
REP.:		CHRISTIAN DEGRACE	
CUSTOMER:		---	
DATE (DM/Y):		13-02-2004	



Curved Track installed through door header and into washroom. Due to the difference in height from the ceiling in both areas the tracking was dropped to accommodate the lower of the two ceiling heights



**Kwiktrak
Installations**

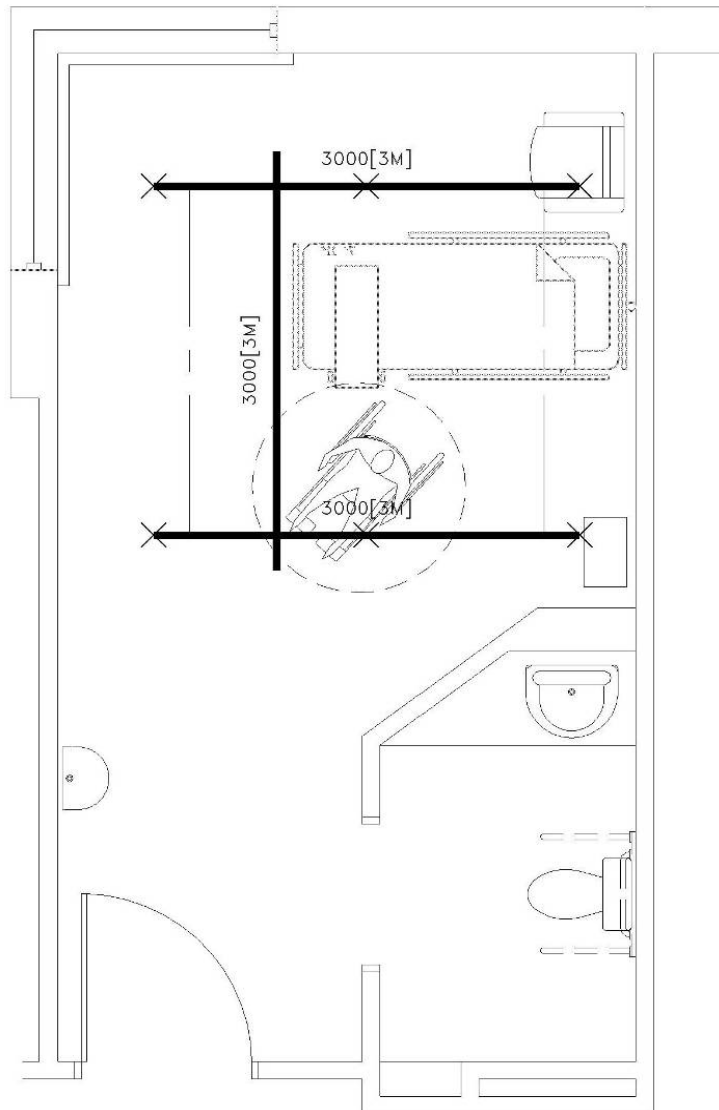

LIST OF MATERIAL		
QTY	PIECE	DESCRIPTION
3	200.12810	3 METERS KWIKTRAK
6	700.11000	END STOPPER
6	200.11060	END CAP
6	700.11100	TRACK BRACKET
1	700.11290	X-Y TROLLEY KIT
1	700.15510	X-Y INTERFACE CHARGER
1	9100201	V4 - 2 WAYS


Private Room XY System:


Configuration Advantages:

- Full floor coverage underneath side rails
- Multiple transfer points
- Use of lift for repositioning and turning in bed
- Installed price guide £2250-2500

Approved: _____

QTY: 1

		DISTRIBUTOR:	
		BHM MEDICAL INC.	
REFERENCE #:		REP.:	
B040213		7	
REVISION:		CUSTOMER:	
---		---	
DRAWN BY:		DATE (D/M/Y):	
F. GUERTIN		13-02-2004	
SCALE:			
1=50			

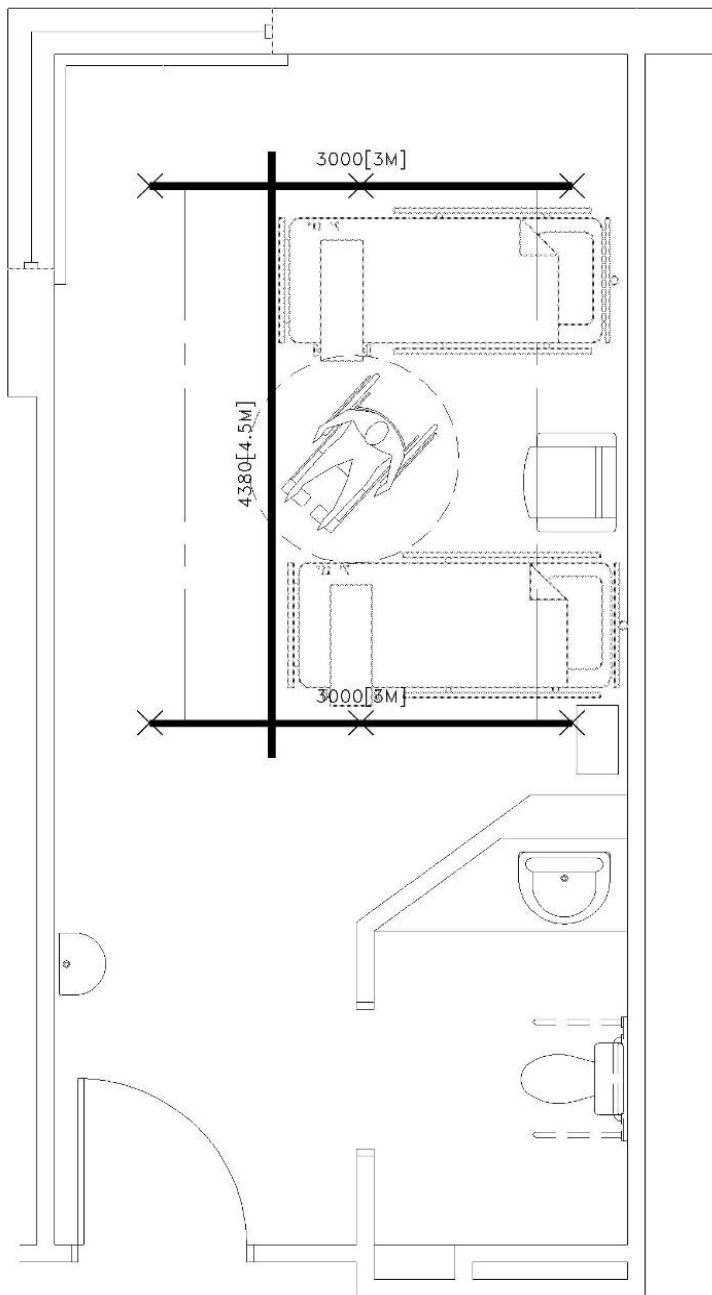

LIST OF MATERIAL		
QTY	PIECE	DESCRIPTION
2	200.12810	3 METERS KWIKTRAK
1	200.12870	4.5 METERS KWIKTRAK
6	700.11000	END STOPPER
6	200.11060	END CAP
6	700.11100	TRACK BRACKET
1	700.11290	X-Y TROLLEY KIT
1	700.15510	X-Y INTERFACE CHARGER
1	9100201	V4 - 2 WAYS


Semi-Private Room XY System:


Configuration Advantages:

- One lift covers both beds
- Full floor coverage underneath side rails
- Multiple transfer points
- Use of lift for repositioning and turning in bed
- Installed price guide £2250-2600

Approved: _____

QTY: 1

 BHM MEDICAL INC.		DISTRIBUTOR: BHM MEDICAL INC.	
REFERENCE #: B040213		6	
REVISION: ---		REP.: CHRISTIAN DEGRACE	
DRAWN BY: F. GUERTIN		SCALE: 1=50	
		DATE (D/M/Y): 13-02-2004	

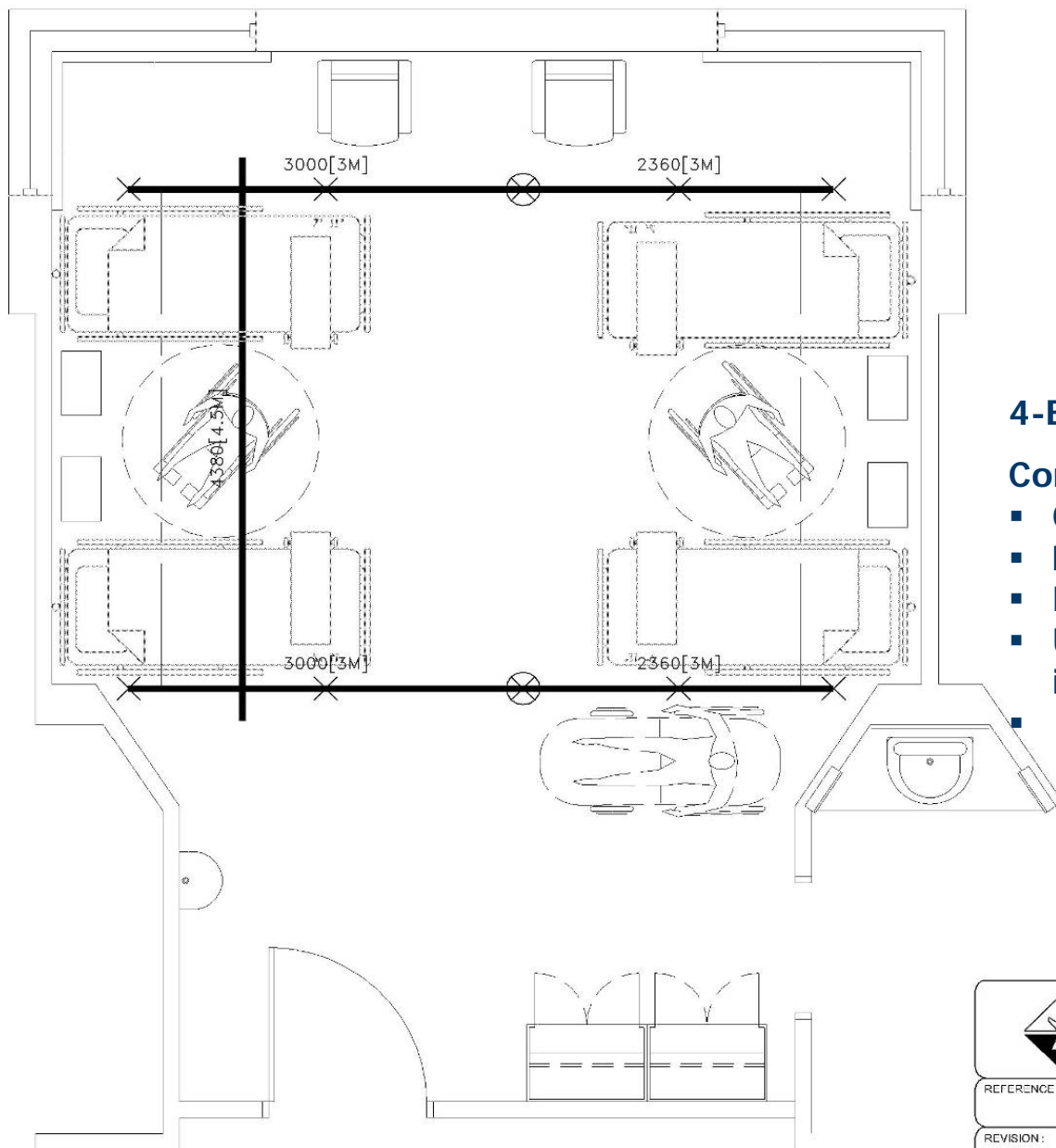

LIST OF MATERIAL		
QTY	PIECE	DESCRIPTION
4	200.12810	3 METERS KWIKTRAK
1	200.12870	4.5 METERS KWIKTRAK
6	700.11000	END STOPPER
6	200.11060	END CAP
10	700.11100	TRACK BRACKET
4	200.11105	CAM FOR JOINT
4	000.00810	SPRING PIN
1	700.11290	X-Y TROLLEY KIT
1	700.15510	X-Y INTERFACE CHARGER
1	9100201	V4 - 2 WAYS


4-Bed Ward XY System:

Configuration Advantages:

- One lift covers all beds
- Multiple transfer points
- Full floor coverage underneath side rails
- Use of lift for repositioning and turning in bed
- Installed price guide £3000-3250

Approved: _____

QTY: 1


 BHM MEDICAL INC.		DISTRIBUTOR: BHM MEDICAL INC.	
REFERENCE #: B040213		5	
REVISION: ---		CUSTOMER: ---	
DRAWN BY: F. GUERTIN		SCALE: 1=50	
		DATE (D/M/Y): 13-02-2004	



XY System installed into a
Private Room




**Kwiktrak
Installations**

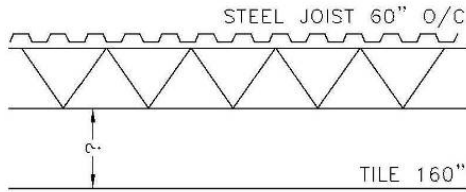

16	200.11105	CAM FOR JOINT
16	000.00810	SPRING PIN
36	700.11350	LATERAL BRACE
4	700.11600	TURNTABLE
4	9100001	V4 - 4 WAYS

YOU SHOULD FORESEE ELECTRICAL
OUTLET FOR CHARGER.

LIST OF MATERIAL		
QTY	PIECE	DESCRIPTION
8	200.12800	2 METERS KWIKTRAK
8	200.12810	3 METERS KWIKTRAK
4	700.12840	45° CURVE KWIKTRAK
12	700.11005	END STOPPER
12	200.11060	END CAP
52	700.11100	TRACK BRACKET

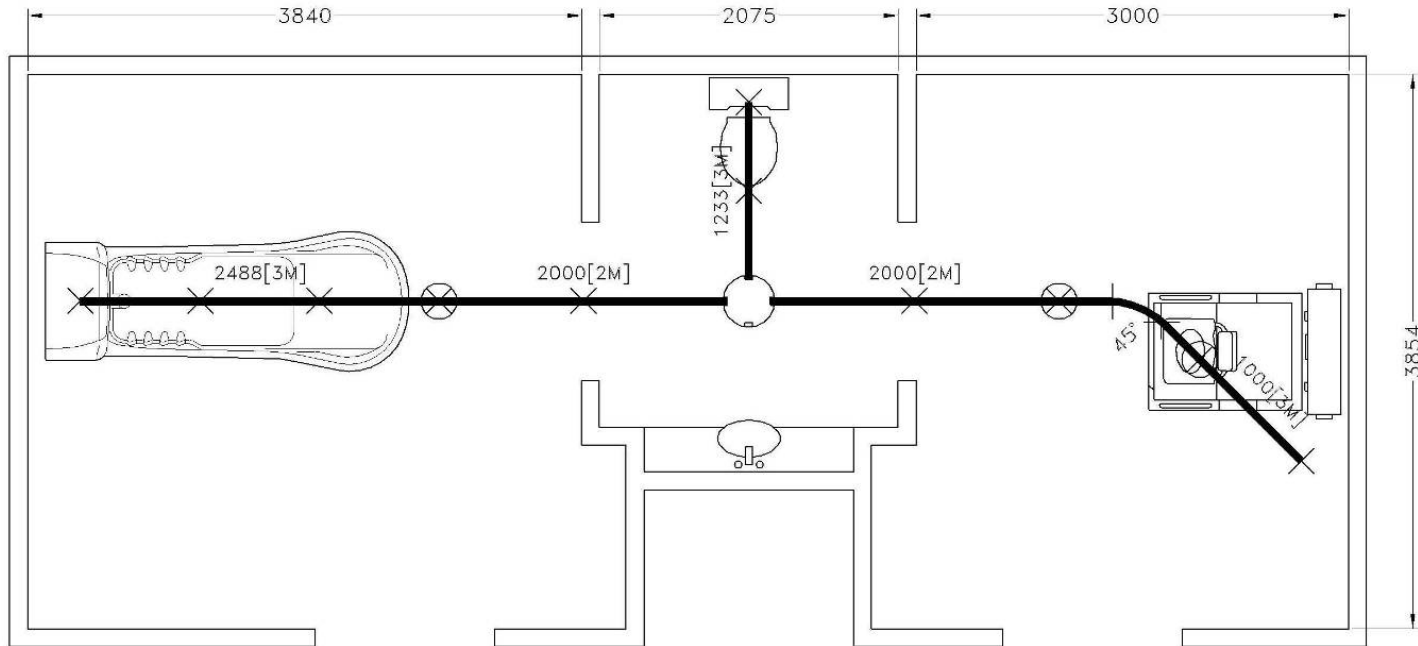

Approved: _____

Segregated Bathroom:


Configuration Advantages:

- With use of turntables one lift covers bathtub, toilet and shower cabinet

APPROXIMATE TRACK LENGTH
(DO NOT CUT TRACK BEFORE INSTALLATION)

This layout has been prepared with uncompleted inputs.
Verifications must be done on site.

PLEASE CONFIRM
DIMENSIONS


		DISTRIBUTOR:	
		REP.:	
REFERENCE #:	1/1	CUSTOMER:	
REVISION:		DATE (D/M/Y):	
DRAWN BY:	SCALE: 1=50		



A turntable allows you to join many tracks running in different directions together in one junction.



**Kwiktrak
Installations**

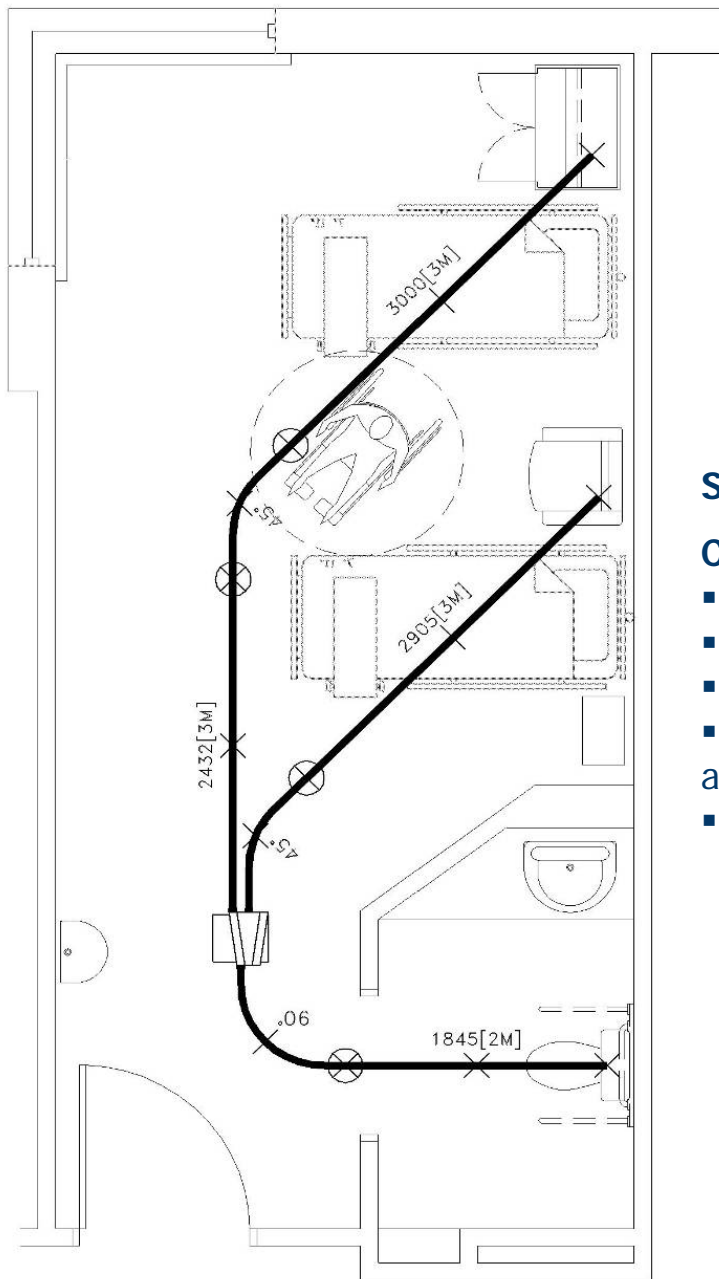

LIST OF MATERIAL		
QTY	PIECE	DESCRIPTION
2	200.12800	2 METERS KWIKTRAK
2	200.12810	3 METERS KWIKTRAK
2	700.12840	45° CURVE KWIKTRAK
1	700.12860	90° CURVE KWIKTRAK
3	700.11000	END STOPPER
3	200.11060	END CAP
14	700.11100	TRACK BRACKET
8	200.11105	CAM FOR JOINT
8	000.00810	SPRING PIN
1	700.11700	EXCHANGER
1	9100001	V4 - 4 WAYS


Semi-Private Room Curved Track with Exchanger:

Configuration Advantages:

- One lift covers both beds
- Multiple transfer points
- Access to washroom with lift
- 4-way lift allows powered movement left and right along rail
- Use of lift for repositioning and turning in bed

Approved: _____

QTY: 1


 BHM MEDICAL INC.		DISTRIBUTOR:	
		BHM MEDICAL INC.	
REFERENCE #:		REP.:	
B040213		CHRISTIAN DEGRACE	
REVISION:		CUSTOMER:	
---		---	
DRAWN BY:		SCALE:	
F. GUERTIN		1=50	
		DATE (D/M/Y):	
		13-02-2004	



An exchanger allows you relocate the lift from one track to another.



**Kwiktrak
Installations**

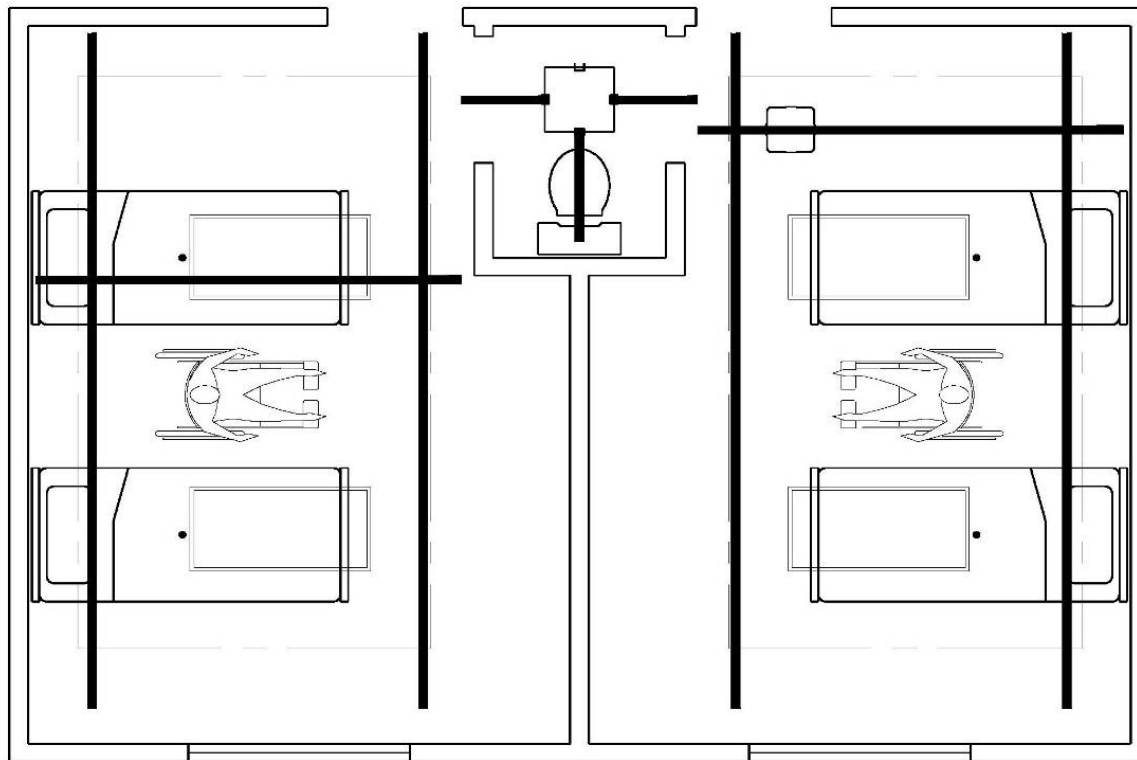

XY System to Gate

Configuration

Advantages:

The addition of a gate to the XY System allows us to transfer the lift from the traverse rail to another set of tracks. This allows us to access bathrooms and other rooms with the same lift.



BHM
MEDICAL INC.


TRACK LAYOUT



A Gate allows you to travel the lift from an XY System to either straight or curved tracking to gain access into washrooms or other rooms.



**Kwiktrak
Installations**

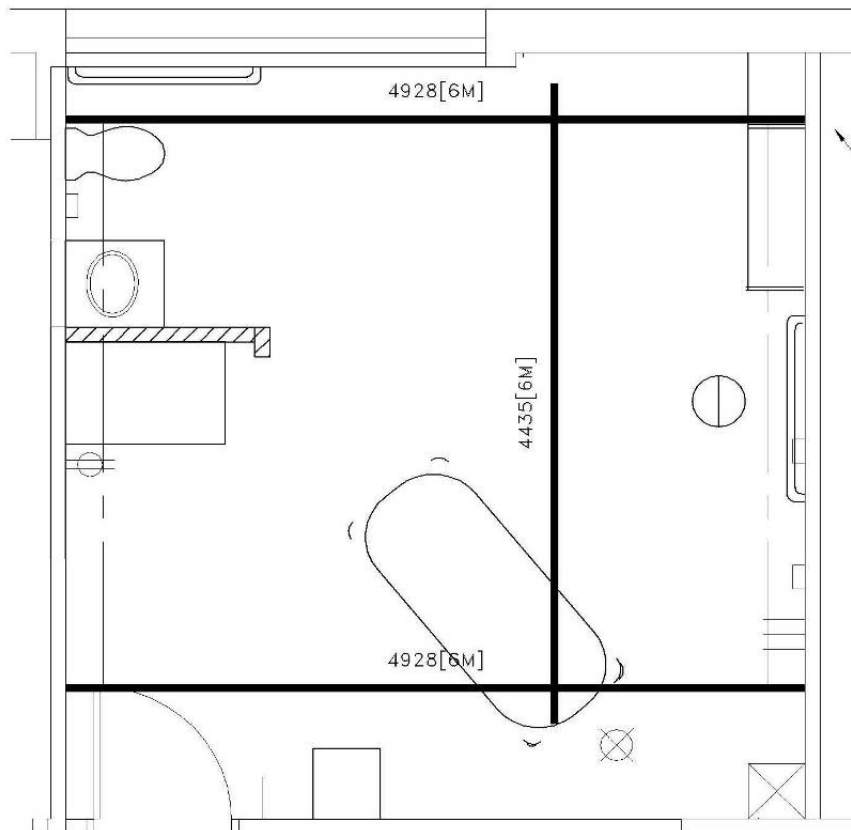

LIST OF MATERIAL		
QTY	PIECE	DESCRIPTION
3	200.12880	DOUBLE 6 METERS KWIKTRAK
1	700.11290	X-Y TROLLEY KIT
1	700.15510	X-Y INTERFACE CHARGER
1	700.11020	WALL BRACKET
3	700.12900	WALL POST KIT
1	9100201	V4 - 2 WAYS


YOU SHOULD FORESEE ELECTRICAL
OUTLET FOR CHARGER.

Wall to Wall Installation:

Configuration Advantages:

- When there isn't an opportunity to install into the ceiling, it is still possible to install tracking with wall brackets and wall posts
- Installation price add £200 to standard install equivalent

Approved: _____

QTY: 1			
 BHM MEDICAL INC.		DISTRIBUTOR:	
REFERENCE #:		1/1	
REVISION:		REP.:	
		CUSTOMER:	
DRAWN BY: F. GUERTIN		SCALE: 1=50	DATE (D/M/Y):

APPROXIMATE TRACK LENGTH
(DO NOT CUT TRACK BEFORE INSTALLATION)

This layout has been prepared with uncompleted inputs.
Verifications must be done on site.

**PLEASE CONFIRM
DIMENSIONS**



Wall to Wall install using wall posts. This is necessary when installation into the ceiling is not possible.














**Kwiktrak
Installations**

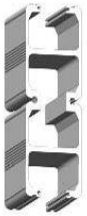

A pivoting curtain rod allows you to drop the privacy curtain underneath the tracking in order to eliminate any interference with the lift.

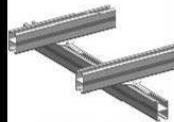


**Kwiktrak
Accessories**

	 Track height 90 mm	 Track height 140 mm	 Track height 180 mm	 Curved sections (with 585 radius)	 Turntable	 Exchanger	 Gate	 Traverse system (track not included)
 MAXI V3 40	Available in 2 m and 3 m lengths. 2 m #200.12800 3 m #200.12810 Max. unsupported span: 2.5 m	Available in 4.5 m lengths. #200.12870 Max. unsupported span: 4.5 m	Available in 6 m lengths. #200.12880 Max. unsupported span: 6 m	Available in 45 and 90 degree. #700.12840 (45°) #700.12860 (90°)	Pull cord with external power supply included #700.11651=UL #700.11652=JP #700.11653=CE #700.11654=AU #700.11655=UK	UNDER DEVELOPMENT	UNDER DEVELOPMENT	2 Trolley Assembly #700.11290
 MAXI V4 100	Available in 2 m and 3 m lengths. 2 m #200.12800 3 m #200.12810 Max. unsupported span: 2.15 m	Available in 4.5 m lengths. #200.12870 Max. unsupported span: 4.15 m	Available in 6 m lengths. #200.12880 Max. unsupported span: 5.95 m	Available in 45 and 90 degree. #700.12840 (45°) #700.12860 (90°)	Powered by the lift cassette #700.11600	Powered by the lift cassette #700.11700	Powered by the lift cassette #700.11500	2 Trolley Assembly #700.11290 Charger interface not included with the assembly #700.15510
 MAXI SKY 1000	X NOT PERMITTED	Available in 4.5 m lengths. #200.12870 Max. unsupported span: 2 m	Available in 6 m lengths. #200.12880 Max. unsupported span: 3.5 m	X NOT COMPATIBLE	X NOT AVAILABLE	X NOT AVAILABLE	X NOT AVAILABLE	X NOT COMPATIBLE

MAX. UNSUPPORTED SPAN: The horizontal distance between two support brackets.



Kwiktrak Matrix

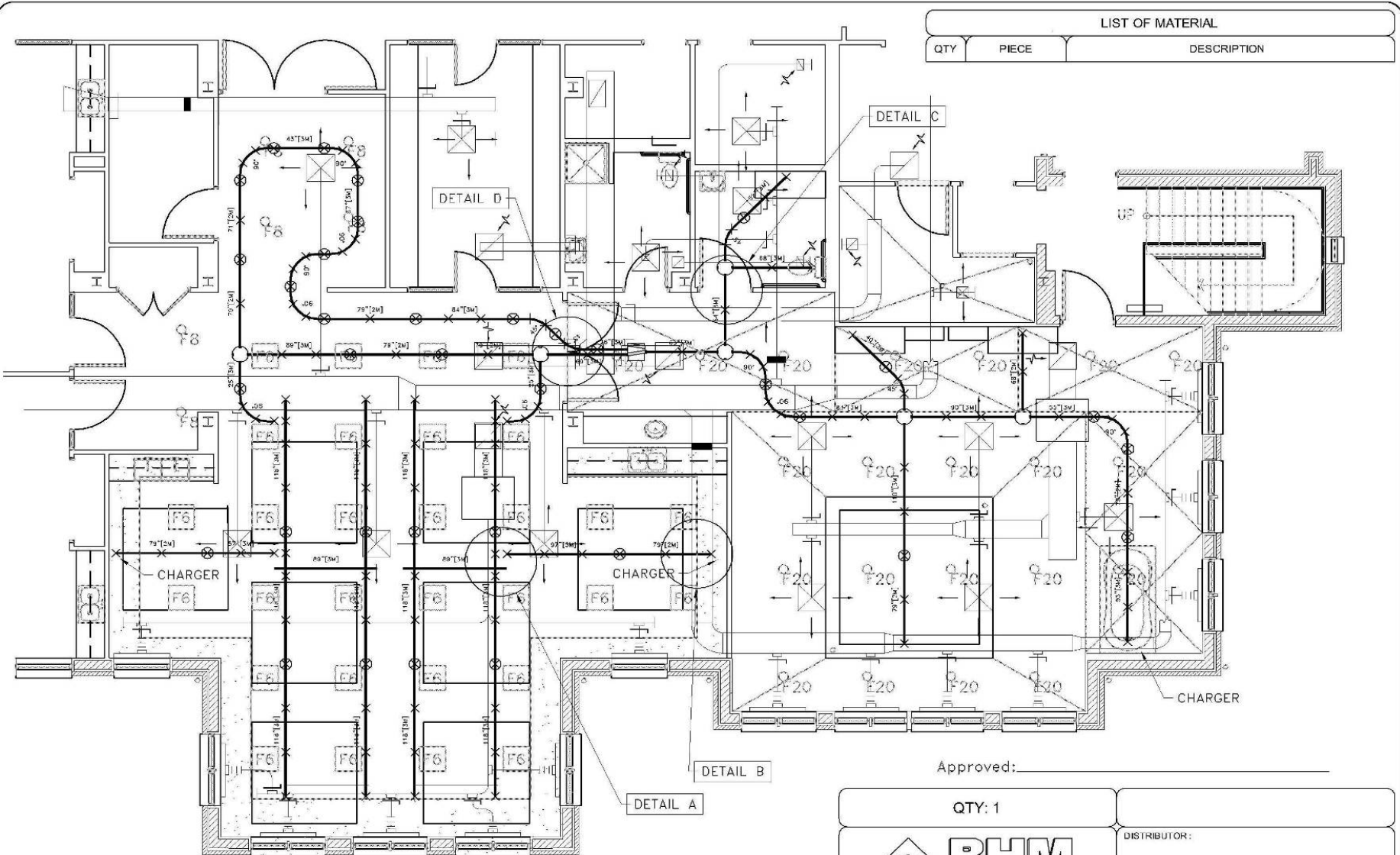

LIST OF MATERIAL

QTY	PIECE	DESCRIPTION
-----	-------	-------------

Anything and everything is possible! We have never failed to find a solution!

APPROXIMATE TRACK LENGTH
(DO NOT CUT TRACK BEFORE INSTALLATION)

This layout has been prepared with uncompleted inputs.
Verifications must be done on site.

PLEASE CONFIRM
DIMENSIONS

Approved: _____

QTY: 1



BHM
MEDICAL INC.

DISTRIBUTOR:

REFERENCE #:

1.1

REP.:

REVISION:

CUSTOMER:

DRAWN BY:

SCALE:

DATE (MM/YY):

1=1/4



Where do we use Ceiling Lifts...



Hospitals, LTC, Rehab Centers, Homecare, Riding Schools: the possibilities are endless



Where ?

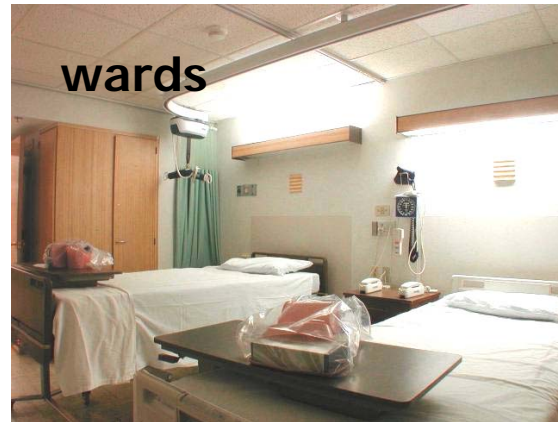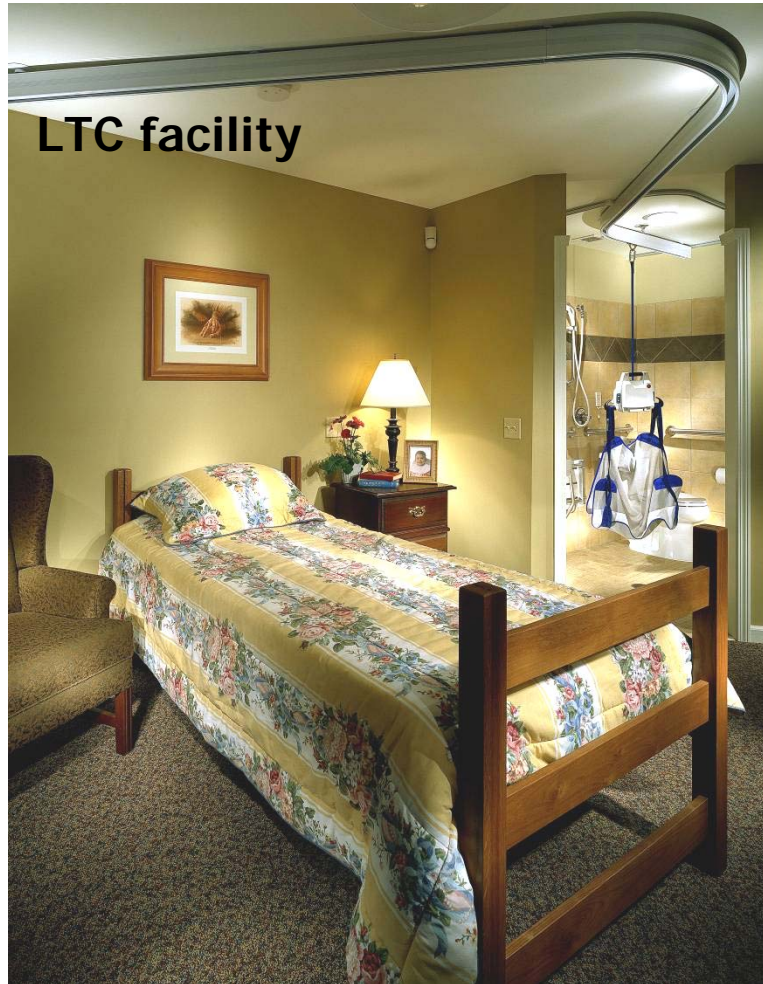

Where ?



Where ?



BHM

Because Health Matters

What configuration is best

- Straight track
- J track
- U track
- X-Y gantry

.... the simpler the better !

